

FRC Postec® Plus

1

Preparation of the post space | Präparation des Stiftbettes

Gates Glidden /
Peeso





opt.
1000-5000 rpm


Size 0: Ø max. 0.5mm

Size 1: Ø max. 0.7mm

Size 3: Ø max. 0.9mm

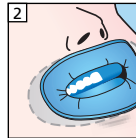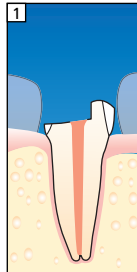

Isolation
Absolute
Trockenlegung


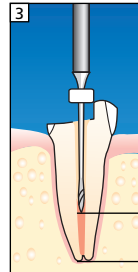


min. 4 mm


FRC Postec Plus
Reamer



opt.
1000-5000 rpm

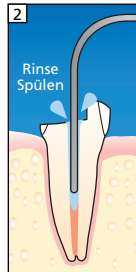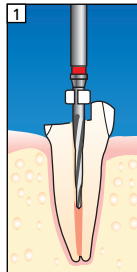


Rinse
Spülen


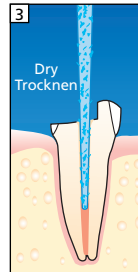


Dry
Trocknen

FRC Postec Plus

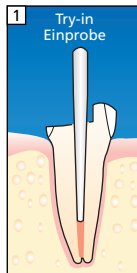


Try-in
Einprobe

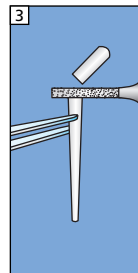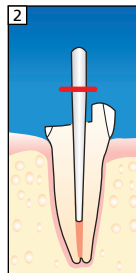


2 Conditioning of the post | Konditionierung des Stiftes

Total Etch

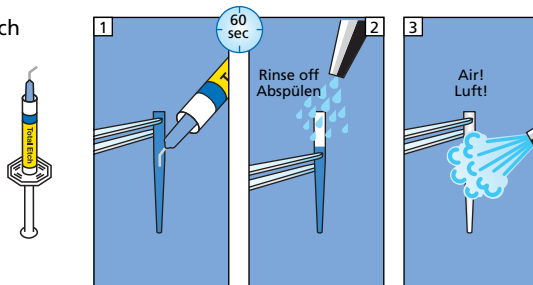

Monobond Plus

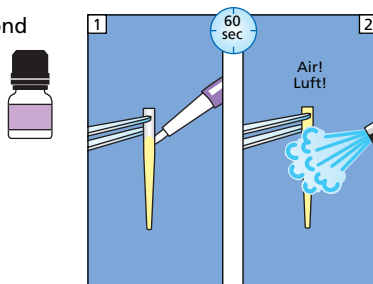

3 Luting of post and Core build-up | Stiftbefestigung und Stumpfaufbau

Please observe the Instructions for Use of the luting composite and the core build-up material.

Bitte beachten Sie die Gebrauchsinformation des Befestigungscomposites und des Stumpfaufbaumaterials.



Please read Instructions for Use
Bitte lesen Sie die Gebrauchsinformation

CE 0123